


INTENSE PULSED LIGHT TECHNOLOGY FOR PHOTOREJUVENATION AND PROGRESSIVE PERMANENT EPILATION



Intuitive touch screen with preset protocols and customizable programs.



Personalization of treatment.

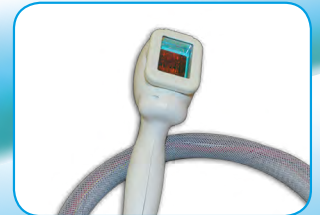

Extractable hand piece with spot of 4,5x3,5 cm.



Also areas of treatment more delicate due to the reduction gear



Filters for:

- Epilation dark skin
- Epilation clear skin
- Photorejuvenation
- Acne
- Vascular

ATHON2

INTENSE PULSED LIGHT



Very comfortable treatment for the operator thanks to the pedal.

ACTIONS:

- PHOTOREJUVENATION • VASCULAR • COUPEROSE • ACNE
- SKIN SPOT • PROGRESSIVE PERMANENT HAIR REMOVAL

ATHON2

INTENSE PULSED LIGHT

THE TECHNOLOGY:

ITS Group research and development department created ATHON 2, the no coherent pulsed light system which evolve the concept of epilation. The technical specifications allow to do treatments of epilation, photorejuvenation, vascular, acne and skin spot in a safety and effective way.

ATHON 2 is the intense pulsed light device most required by the medical centers for his possibility to personalize the treatments by setting intensity and pulse time, and so you can use the correct fluency in joule/cm² for any case.

ATHON 2 allows to treat big areas in a short time thanks to the spot of 45x35mm. The wide range of filters supplied, which are different in their double coating allows to treat any kind of hair and skin.

ATHON 2 is able to act on the microcapillaries that bring nutrients to the germinative papilla thus causing atrophy of hair bulbs, for this is able to treat even white hairs and dark skins.

Visible and lasting results already after 3-4 sessions.

CONTROINDICAZION:



- Epilepsy
- Menstrual cycle in progress
- Last sun exposure 45 days
- Taking Cortisone
- Taking antibiotics
- Taking vitamin A e K
- Presence of perfumes or potentially irritating substances
- Pregnancy
- Pace maker
- Tan
- Skin lesion

INFORMATION:



- Number of sessions : 6/12
- Session interval: once or twice per month
- Session duration : 5/50 minutes depending on the area

TECHNICAL SPECIFICATIONS:

Model:	ATHON2
Class:	1
Supply:	220-240V 50-60Hz
Power consumption:	400W
Fluency:	max 70J/cm ²
Spot size:	35x45 mm (1575 mm ²)
Charging time:	~2 sec
Pulse duration (ms):	4-50 ms
Standard wavelength:	590nm-600nm – 650nm
Classification IEC 60601-1:	Class I type B
Classification 93/42/CEE:	II a
Display:	Touch Screen
Weight:	17,5 kg
Size:	L 450 x P 500x H 250 mm
Conformity and Regulations:	Temperature environment: 15°C - 35°C, moisture 30% - 60%, pressure 86kPa -106kPa
Storage and transport:	Temperature environment: 5°C- 45°C, moisture 30% - 70%, pressure 86kPa -106kPa



PERMANENT EPILATION
CLEAR SKIN



PERMANENT EPILATION
DARK SKIN



SKIN SPOT



PHOTOREJUVENATION



ACNE



VASCULAR

