

## OXIGEN DERMO INFUSION FOR SKIN AND SCALP TREATMENT



Specific treatment for scalp:  
dandruff, purifying, sensitive  
skin, mineralizing and fall



The only non-invasive treatment  
that allows it to operate on the up-  
per eyelid and the eyeline.



Lighting and oxygenation of  
the skin.



Applications to 360°, even on  
large areas.



Equipped with a special  
airbrush with adjustable flow.



Special MICROCAMERA to analyse  
the skin and the sculp with the  
possibility to take pictures.

## OXYPROAGE PLUS OXYGEN DERMO INFUSION

### ACTIONS:

#### FACE TREATMENTS:

- DARK CIRCLES • EYE CONTOUR
- FACE LIFTING • ACNE • WRINKLES • SPOTS

#### SCALP TREATMENT:

- ANTI HAIR LOSS • REMINALIZER
- SENSITIVE SKIN • PURIFYING
- ANTI DANDRUFF



Nono technological serums  
specifically created for for the  
transdemc vehiculation to treat  
every blemish of the skin.

MADE IN ITALY



# OXYPROAGE PLUS

## OXYGEN DERMO INFUSION

### THE TECHNOLOGY:

Oxygen is precious and essential to the life of each of the cells of our body. A natural and perfectly biocompatible component that favors metabolic exchanges within the tissues is now available in a pure state to be applied in the aesthetic treatments at the highest level.

Designed and tested by the research and development department of ITS Group, OXYPROAGE PLUS is the new and sophisticated technology that combines the Oxygen Dermo Infusion with the nanostructured actives. The pure oxygen filtered and delivered to the skin in continuous propulsion of 3 liters per min has the property to bind with the nanostructured active principles and convey them in depth, below the stratum corneum, to carry out their function in the target tissue adequately prepared. The synergistic actions of the oxygen dermo infusion and nanostructured serum give a unique result, immediate and long lasting. This advanced combination allows to get quick and brilliant results on the tissue: it stimulates tissues oxygenation and microcirculation and the skin will be immediately more bright, smoothed, superficial wrinkles plumped, uniform, less visible pores, attenuation of dark circles and the skin will be visibly cared and nourished.

We will also have an anti-inflammatory action on the acne and a lightening effect on the hyperchromatic spots. The Oxygen Dermo Infusion stimulates tissues oxygenation and microcirculation, revitalizing the cutis and enhancing hair shine and vitality. The special scalp serums can solve the main problems of the cutis (oily and heavy hair, seborrhea, dandruff, thin and devitalized hair, sensitive skin, hair loss and thin hair). The device is also equipped with a light and handy digital microcamera, which is ideal for dermatological and trichological analysis before, during and after treatment. The microcamera is managed by an easy and intuitive software already installed on the device, which can take and save pictures to record client's progress. The microscope is also equipped with a dock complete with a USB door for charging batteries via PC.

### CONTRAINDICATIONS:




- Pregnancy
- Pacemaker

### INFORMATION:



- Number of sessions: 4/12
- Session interval : once or twice a week
- Session duration : 30/40 minutes

### OPTIONAL

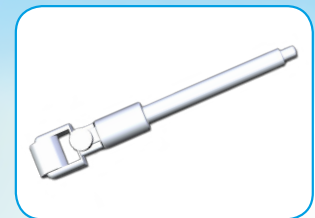

Applicator roller with needles to enhance the penetration of actives.



### SKIN SERUMS

1. **AGE REPAIR SERUM:** *Anti-Wrinkle, Elasticity and Hydrating Action*
2. **CRYSTALL FILLER SERUM:** *Filler Action*
3. **DEFENCE SYSTEM SERUM:** *Lenitive, Anti-Inflammatory And Restructuring Action*
4. **PERFECT SKIN SERUM:** *Astringent, Dermopurifying Action*
5. **ACTIVE LIGHT SERUM:** *Lightening Action*
6. **LIPOACT SERUM:** *Lipolytic And Draining Action*

### SCALP SERUMS

1. **APPLE STEM SERUM:** *Anti Hair Loss*
2. **SUNFLOWER SERUM:** *Protective and Remineralizer*
3. **MARIGOLD SERUM:** *Sensitive and Irritated Skin*
4. **LARCH SERUM:** *Purifying Serum for Sebaceous Hyperplasia*
5. **KELP SERUM:** *Anti Dandruff*
6. **OXYGEN BOOSTER:** *Regenerating and Revitalizing*

### TECHNICAL SPECIFICATIONS:

Model:	Oxy Pro Age
Supply:	110 – 240 V; 50-60 Hz
Power consumption:	400W
Oxygen flow:	1 – 5 l/min
Output pressure:	50–80 kPa
Concentrazione ossigeno:	93±3%
Classification IEC 60601-1	Class I type B
Weight:	21 kg
Size:	L380 x P340 x H620 mm
Compliance and Standards:	Temperature environment: 15°C - 35°C, moisture 30% - 60%, pressure 86kPa -106kPa
Storage and transport:	Temperature environment: 5°C- 45°C, moisture 30% - 70%, pressure 86kPa -106kPa

